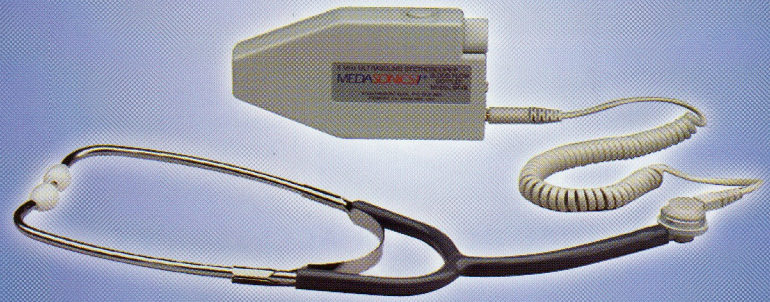


# Model BF4B/5A Handheld Doppler



- **Leader**

*The preferred blood flow Dopplers for assessing peripheral vascular disease*

- **Reliable**

*Durable design performs even in the most rugged environments*

- **Classic**

*The Gold Standard for outstanding hand held Doppler performance, featuring narrow beam technology for increased sensitivity*

- **Adaptable**

*Special Model BF5A available with precision probe tip to easily identify and assess blood flow in small superficial vessels*

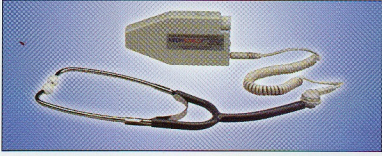
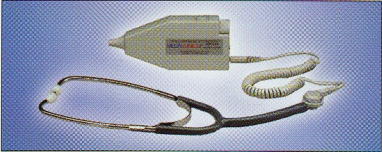
**Model BF4B/BF5A Blood Flow Doppler**



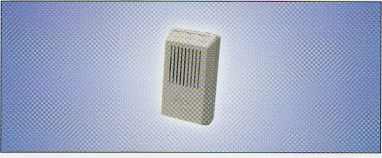
# Model BF4B & BF5A

# Blood Flow Doppler Ultrasound Stethoscope®

## ORDER INFORMATION

	Description	Order Number	Frequency	Contents - Applications / Specifications
	Model BF4B General Blood Flow Doppler	101-0154-010	5.3 MHz	Includes: 5 MHz Doppler, Stethoscope headset, gel, battery, soft carrying case. Applications: Locate vessels and assess vessel patency, calculate ankle-brachial index, measure systolic blood pressure.
	Model BF5A Superficial Blood Flow Doppler	101-0137-010	8 MHz	Includes: 8 MHz Doppler, Stethoscope headset, gel, battery, soft carrying case. Applications: Locate superficial vessels; skin flap transfer, routine arterial cannulations, vascular grafting and microvascular surgery.

## OPTIONAL DOPPLER ACCESSORIES

	Model SA2 Pocket Speaker	101-0052-030	n/a	Includes: Speaker, cord to connect to Doppler. Specifications Power Source: one 9V alkaline battery Size: 4.2 x 2.5 x 1.1 in (10.8 x 6.3 x 2.7 cm) Weight: 6 oz (170 gm)
	Model SA3 Tabletop Speaker	101-0009-010	n/a	Includes: Speaker, cord to connect to Doppler. Specifications Power Source: four 9V alkaline batteries Size: 7.9 x 6.7 x 5.1 in (20 x 17 x 13 cm) Weight: 3.75 lb (1.7 kg)
	Model SA6 Tabletop Speaker	101-0240-010	n/a	Includes: Speaker, cord to connect to Doppler. Specifications Power Source: four C-cell alkaline batteries Size: 5.5 x 3.5 x 6.5 in (15 x 8.9 x 16.5 cm) Weight: 2.25 lb (1 kg)

**Additional Accessories:** Replacement cords, stethoscopes and repair kits are also available. Contact your Sales Representative for part numbers and pricing.

## Model BF4B & BF5A SPECIFICATIONS

## DISTRIBUTED BY

<b>Ultrasonic Frequency</b>	BF4B: 5.3 MHz BF5A: 8 MHz
<b>Outputs</b>	Two 2.5 mm jacks to drive one or two headsets, one headset and a speaker or recording device.
<b>Controls</b>	Volume control knob and push button ON/OFF
<b>Power Source</b>	One standard 9V alkaline battery
<b>Weight</b>	4.7 oz. (133 gm)
<b>Dimensions</b>	4.95 x 2.5 x 1.0 in (12.6 x 6.4 x 2.5 cm)
<b>Construction</b>	Rugged ABS plastic case, integral transducers, solid state circuitry
<b>Service:</b>	MedaSonics, Inc. maintains a full service facility which can promptly repair products returned to the factory.
<b>Warranty:</b>	MedaSonics Dopplers are covered by a 5-year warranty on parts; 1-year warranty on labor. Call Customer Service for details.

## MANUFACTURED BY

**MEDASONICS//®**  
 38875 Cherry Street Newark, CA 94560 USA  
 510/494-1097 or 800/227-8076  
 fax 510/494-0448  
 email: info@MedaSonics.com