



ACCUTEMP

**ACCURACY....THE HEART OF WHAT'S BEEN MISSING
IN TEMPERATURE MEASUREMENT.**



NON IRRITATING



EASE OF USE



SAFER ADHESIVE



LOW PROFILE



CONVENIENT DISPENSER



ECONOMICAL

ACCUTEMP

Clinical investigations of infants in incubators and radiant warmers clearly show the importance of accurate temperature maintenance in order to properly manage apnea, metabolic rate and insensible water loss. Our new ACCUTEMP helps achieve this goal by solving the problems associated with skin probe placement! Here's how:


The insulating foam layer keeps the ambient air from affecting the temperature sensor's determination of patient skin temperature.

♥ Gold reflective mylar foil further shields the probe from the radiant heat effects of overhead warmers. With ACCUTEMP you're assured the probe is being warmed by the patient...not the warmer.

♥ The temperature sensor is held securely in place to optimize heat conduction from the skin to the probe.

♥ ACCUTEMP's "easy grip" adhesive simplifies removal without leaving a messy adhesive residue on the patient or the probe.

♥ Two convenient sizes are offered; Choose the one which best meets your clinical needs.

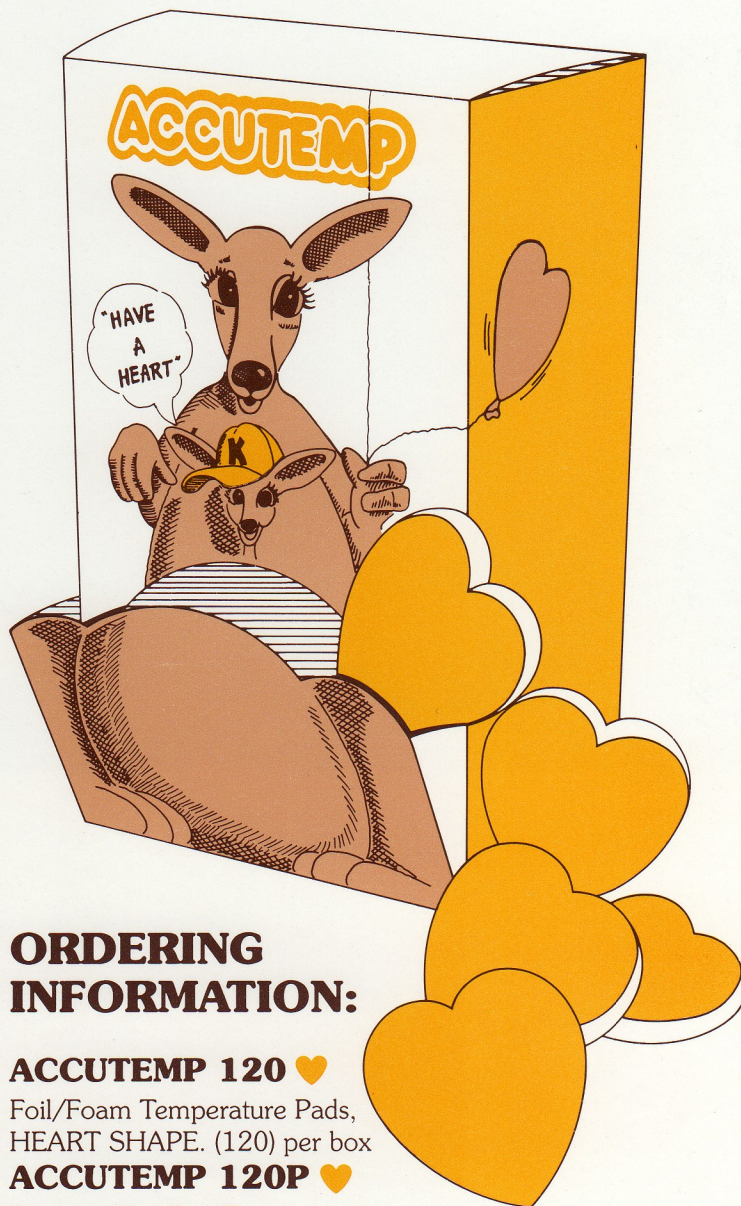

ORDERING INFORMATION:

ACCUTEMP 120 ♥
Foil/Foam Temperature Pads,
HEART SHAPE. (120) per box

ACCUTEMP 120P ♥
Premie Size Foil/Foam
Temperature Pads,
HEART SHAPE. (120) per box

DISTRIBUTED BY: